

Wart Treatment Liquid

Carefully read the instructions for use prior to using Dr. Yglo Wart Treatment Liquid, even if you used the product before. The information may have been amended according to new research and findings in the time since you opened your last package. Keep these instructions for use in case you need to consult them again. For adults and children (4 years and over).

What is Dr. Yglo Wart Treatment Liquid used for?

Dr. Yglo Wart Treatment Liquid is a medical device for the treatment of common warts. Dr. Yglo Wart Treatment Liquid contains monochloroacetic acid, a substance which has a keratolytic effect on the skin. The wart will dissolve, layer by layer. Monochloroacetic acid ensures that the cells containing the wart virus are destroyed.


The benefit of using Dr. Yglo Wart Treatment Liquid is that an application directly onto the wart does not hurt, which means it can also be used for children, when applied correctly by an adult.

Dr. Yglo Wart Treatment Liquid can only be used on common warts:




- Common warts primarily occur on hands and feet
- Same color as the skin
- Circumference varies from less than a millimetre to more than a centimetre
- Several warts usually occur simultaneously
- Contamination often occurs through direct contact, but also through towels or infected surfaces
- Scratching and nail biting increase the risk of contamination

Do not use Dr. Yglo Wart Treatment Liquid on senile warts:

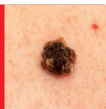

- Senile warts can come in different sizes and colors
- Average circumference is 0.5 to 1 cm
- Often start as a skin-colored, yellow brown, greyish-brown patch. Later on, they become bigger and fatter
- Appear as if they are lying on the surface of the skin and are easy to scratch off, but are actually firmly attached to the skin
- They can occur anywhere on the body, except for the palms of the hands and soles of the feet, and the lips. Most likely to occur on the chest, the back, the face, and the neck

Do not use Dr. Yglo Wart Treatment Liquid on plantar warts/verrucae:

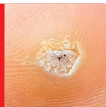

- Plantar warts/verrucae are found on the bottom/soles of your feet
- Usually flat and surrounded by rough, thickened skin
- Resemble calluses
- Can contain black dots or tiny seeds

Instructions for use (Prevent burns):

DO's:

- Application by adults
- Use only on hands and feet. Do not use on other body parts
- Rub Vaseline or zinc ointment onto the surrounding skin to protect the skin
- Place the bottle on a solid underground that is able to resist the caustic fluid
- Dab the wart just once when applying
- Dab the wart only once per week
- Allow the liquid to only come into contact with the wart
- Avoid contact with healthy skin
- Allow the liquid to dry properly after applying
- Leave the wart uncovered
- The result is visible after 1 week
- Keep out of reach of children

DON'Ts:

- Do not apply by individuals under the age of 18;
- Do not apply to warts on any part of the body other than the hands or feet;
- Do not apply on warts that occur completely or partially on the soles of the feet;
- Do not apply to any wart that is wholly or partially under a nail;
- Do not allow the liquid to come into contact with healthy skin;
- Do not apply when your hands are shaking or your fingers are stiff;
- Do not dab on wart more than once when applying;
- Do not cover the wart after applying the liquid;
- Do not apply on birthmarks, water warts, open wounds, hyperpigmentation, senile warts, plantar warts, dark-colored, hairy or weird looking warts;
- Do not apply during pregnancy or lactation.
- Do not repeat treatment without a one-week interval.

⚠ CAUTIONS

Do not use Dr. Yglo Wart Treatment Liquid:

- Under the age of 4 years old
- If you are sensitive to one or more of the ingredients
- On the face (the fluid should never make contact with the eyes or the mucous membranes of the nose or mouth)
- On genital warts or in the genital area
- On birthmarks, molluscum contagiosum, open wounds, moles, naevus, senile (brown) warts, plantar warts, darkly colored, hairy or strange looking warts
- On red, irritated, or infected warts
- Immediately after a surgical removal of warts, freezing of warts (cryotherapy) or a laser treatment. This particularly applies when the skin is red or painful after the treatment. In this case you will have to wait at least 4 weeks, until the redness or inflammation of the skin has fully disappeared. Only then you can start using Dr. Yglo Wart Treatment Liquid again
- During pregnancy or while breastfeeding
- On warts which occur completely or

partially underneath the nail

- Dab wart no more than once and repeat only if necessary after a one-week interval
- Avoid contact with eyes. Contact with eyes may cause blindness
- Do not rub or touch eyes with hands if the product has been applied to them. If product gets into eyes, call medical help right away

Be extra careful when using Dr. Yglo Wart Treatment Liquid:

- In case of doubt concerning the diagnosis (the change in the skin cannot be clearly explained as a wart)
- In case of immune deficiency (e.g. if you are using cytostatic medication)
- In case of diabetes, or other complaints such as vascular diseases, as wounds will heal less easily

In these cases, you may not use Dr. Yglo Wart Treatment Liquid before a doctor has been consulted.

⚠ WARNINGS

- Only intended for external use
- Use with utmost care. The fluid of Dr. Yglo Wart Treatment Liquid is corrosive
- If the healthy skin makes contact with the fluid, please rinse it with water. If, despite these precautions, the skin becomes irritated, it will have to undergo the same treatment as a burn
- Contact your doctor if your complaints do not subside
- Properly close the bottle after use and store it in an upright position
- Keep out of reach of children
- Use by children, the elderly and less-abled must be supervised by an adult
- Avoid contact with clothes, textile, tablecloths or skin
- Stop using if blisters occur
- Dab wart no more than once and repeat only if necessary after a one-week interval

Side effects

Dr. Yglo Wart Treatment Liquid may cause some side effects, although not everyone will experience these side effects. The side effects that could occur when using Dr. Yglo Wart Treatment Liquid are irritated skin, burning sensation, skin redness or

oversensitivity reactions. If any of these side effects become serious, or if you notice any side effects not listed in these instructions for use, please consult your doctor or pharmacist.

Interference with medications

No known interactions with other medication

Contents of Dr. Yglo Wart Treatment Liquid

A bottle filled with with a 2ml mixture of Monochloroacetic acid (50%) and purified water. A custom applicator is used to administer the formulation to the wart.

Storage conditions

Store Dr. Yglo Wart Treatment Liquid at a temperature of 15-25°C in a dry and dark place in the original container and out of reach of children.

The expiry date is stated on the bottom of the packaging. After opening, do not use longer than 4 weeks.

For more information about Dr. Yglo Wart Treatment Liquid visit www.dryglo.com. For all product enquiries, please contact the distributor.

Any serious incident that has occurred in relation to this product should be reported to the manufacturer and to the competent authority of the Member State in which the user is established.

TheOTCLab
Fred. Roeskestraat 115,
1076 EE Amsterdam,
The Netherlands

Company name
Streetname 115,
1076 EE City,
Country

Explanation of symbols

The warning signs and symbols on the packaging and instructions for use are essential to ensure that you use this product safely and correctly and to protect you and others from injury. The symbols are explained in the table below:

Symbol	Explanation
	Symbol for the manufacturer
	Symbol that indicates to consult the instructions for use
	Symbol for the distributor
	Symbol that indicates that caution is necessary when operating the device
	Symbol for the expiration date
	Symbol that indicates that the item is a medical device
	Warning: Corrosive substance
	Not for children aged 3 years and younger
	Symbol that indicates to keep this product out of reach of children
	Symbol that indicates that this product should not be used by pregnant women or while breast-feeding
	Symbol for the Lot number
	Symbol that indicates that the product should be stored between 15°C and 25°C
	Symbol that indicates to keep the product away from sunlight
	Symbol that indicates the Unique Device Identifier

How to use Dr. Yglo Wart Treatment Liquid


(Application only by adults):

Please apply with utmost care.

The solution is caustic. Be extra careful when applying it to parts of your hands or feet that are more difficult to reach. This particularly applies to people with limited mobility or in case of disturbed fine motor functions in the system (trembling hands, stiff fingers).



Apply on
common warts
on hands and
feet only

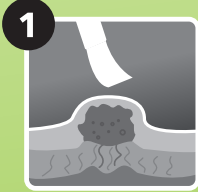

1. Dr. Yglo Wart Treatment Liquid has a child proof cap that can be opened as follows: Place the vial on a solid underground that is able to resist the caustic fluid. Avoid contact with clothes, textile or table-cloths. You can open the vial by pressing down on the cap and turning it counterclockwise simultaneously.
2. Before removing the cap entirely from the vial, apply Vaseline or zinc ointment generously to the healthy skin surrounding the wart.
3. Remove the cap from the vial.
4. Apply a small amount of Dr. Yglo Wart Treatment Liquid to the wart using the applicator (do not allow the fluid to drip). Make sure that you only apply Dr. Yglo Wart Treatment Liquid to the wart and do not dab the wart more than one time.
5. Allow it to dry and do not cover the wart. Carefully close the bottle after use. Dr. Yglo Wart Treatment Liquid should not be applied more than once per week.

The fluid will need a week to dissolve the top layer of the wart, so you will not see an immediate result after application. Depending on the size and thickness of the wart, it may be necessary to repeat the treatment with Dr. Yglo Wart Treatment Liquid after one week. Do not repeat treatment within a one-week interval.

2ml CE 0482



WARNING:
PREVENT BURNS:
CONSULT INSTRUCTIONS
FOR USE BEFORE USING!

